

ScopeTrap

Polyp Trap


Essential device
during all polypectomy
procedures

Tubing with
biconical
connector for the
connection to the
suction device

Tubing for
connection to
the endoscope

Fluorescent
marking of
compartments

Compartments for
resected polyps

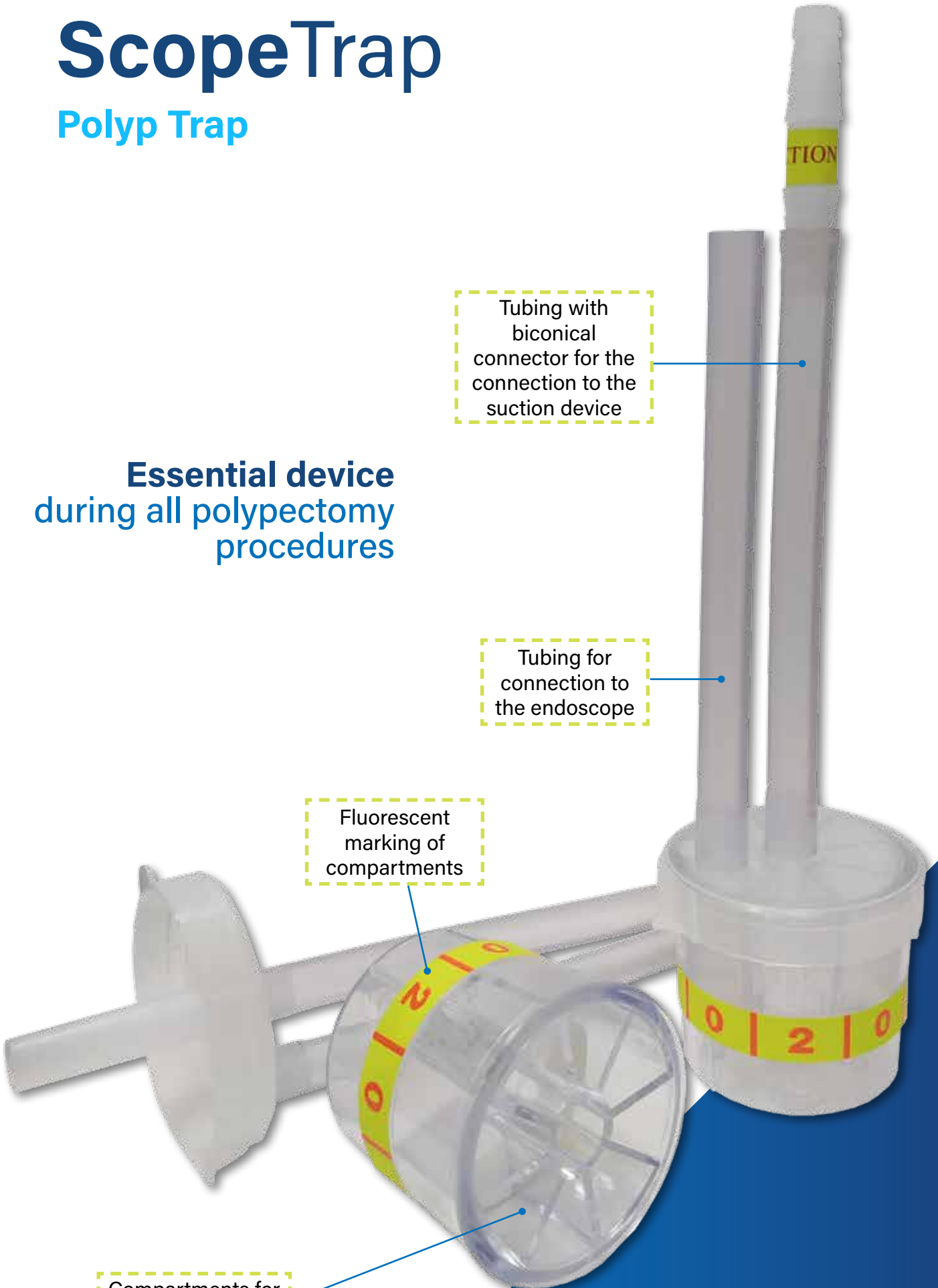

ScopeTrap

Polyp trap

Designed for comfort and security for hospital staff and enabling the filtration, separation and collection of resected polyps by endoscopic means, **ScopeTrap** imposes itself as being an essential device for polypectomy procedures.

ScopeTrap is composed of:

- ▶ A round recipient divided into 8 compartments: 4 (marked from 1 to 4) to chronologically receive the resected polyps, and 4 (marked 0) for the suction.
- ▶ 1 lid containing 2 tubes, one central tube to be connected to the suction device and another tube which is attached to the endoscope.



Reference	Description	Box QTY
10017130	ScopeTrap: Polyp Trap – box of 50 units	50
10017130B	ScopeTrap: Polyp Trap – box of 20 units	20

Model with removable drawers also available, please contact us

Characteristics



- ▶ **Captures retrieved polyps in a closed system**, to maintain the integrity of the samples.



- ▶ **Recipient with 4 numbered compartments**, respecting the chronology of resected polyps.



- ▶ **Fluorescent marking of the compartments** enabling the easy identification of the compartment numbers, even in the half-light of the endoscopy procedure room.



- ▶ **Traceability:** 2 self-adhesive removable labels are included with each polyp trap.

Description: Polyp Trap ScopeTrap. **Indications:** Filtering, separation and collection of resected polyps by endoscopy. **Class:** I. **Certified body:** Auto Certification. **Manufacturer:** Creo Medical S.L. **Use:** Before use, read attentively the instructions in the Instructions for Use or on the product label. **Version:** TRAP_E_3. Rev.03. Date 19-09-22

Creo Medical authorised local distributor

Creo Medical UK Ltd

Chestnut House, N.Minster Business Pk.

YO26 6QR York, United Kingdom

Email: info.uk@creomedical.com

Tel: +44 (0) 1904 786888

Technical changes and mistakes reserved. Images may show non standard options.

