

ScopeValve

Air/Water and Suction channel cleaning adapters

For the pre-cleaning of endoscopes directly after the endoscopic procedure

Air-water Channel Cleaning Adapter

To reduce biofilm formation by removing contaminants in the air and water channels

Biopsy Port Cleaning Adapter

To reduce biofilm formation by removing contaminants in the biopsy port



CREO
MEDICAL

ScopeValve

Air/Water and Suction cleaning adapters

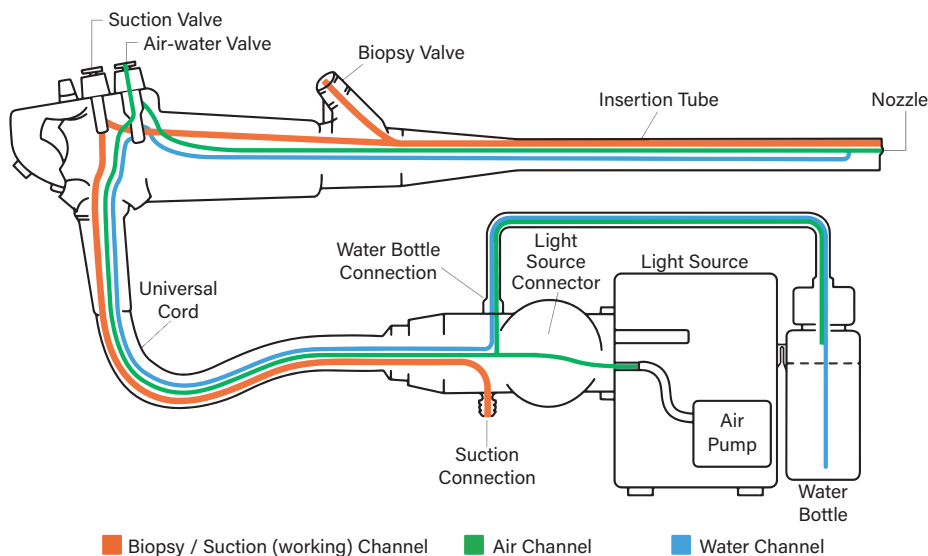
It is strongly recommended to carry out a point-of-use pre-cleaning of each endoscope immediately after the procedure, as the removal of organic substances present after the examination is essential to ensure the effectiveness of the later cleaning steps and hence the final result of disinfection.

ScopeValve air/water channel cleaning adapters are used to rinse the endoscope's air/water channel directly after the patient procedure, preventing biofilm formation and minimizing the risk of cross contamination.





The suction channel cleaning adapter is used to aspirate reprocessing fluids through the endoscope's instrument channel port.

Single use, these devices:

- ▶ Are ready to use,
- ▶ Eliminate manual cleaning and reprocessing of reusable valves and adapters,
- ▶ Reduce the workload of staff responsible for cleaning endoscopes,
- ▶ Improve patient safety with regard to potential contamination risks.



Single use Air/Water and Suction adapters

			
ENDO038	ENDO105	ENDO075	ENDO3875
Reference	Designation		Pkg
ENDO038	Air/Water channel cleaning adapter for Olympus® endoscopes, non sterile		100
ENDO105	Air/Water channel cleaning adapter for Fujifilm® série 700 endoscopes, non sterile		100
ENDO075	Suction channel cleaning adapter for Olympus® endoscopes, non sterile		50
ENDO3875	Suction channel cleaning adapter + air/water channel cleaning adapter for Olympus® endoscopes, sterile		25

Characteristics



▶ **Single Use** for a product which is always effective and a significant reduction in the risk of nosocomial infections.



▶ **Traceability:** 3 detachable self-adhesive labels accompany each device.

Description: Valves and cleaning adapter, ScopeValve. **Indications:** For flushing air/water channels and suction of endoscopes. **Class:** I non-sterile **Notified body:** EC Certificate of Conformity. **Manufacturer:** GA Health Company Limited. **Proper use:** Before use, carefully read the instructions in the package insert or on the device label. **Version :** SVA_E_1. Rev.01. Date 15-03-23

Creo Medical authorised local distributor

Creo Medical UK Ltd

Chestnut House, N.Minster Business Pk.

YO26 6QR York, United Kingdom

Email: info.uk@creomedical.com

Tel: +44 (0) 1904 786888

Technical changes and mistakes reserved. Images may show non standard options.

