

# ScopeXtract

## Triple lumen Biliary Extraction Balloon

Our range of ScopeXtract biliary extraction balloons are designed to remove stones from the bile duct or to facilitate the injection of contrast media.

- ▶ 3 lights: guide wire - injection - balloon inflation
- ▶ 1 balloon = 3 diameters
- ▶ Available in standard or short guide wire model
- ▶ 2.3 mm diameter, tapered to 1.7 mm distally
- ▶ Platoid shaped balloon


**Designed to remove stones from the bile duct or to facilitate the injection of contrast media.**



# ScopeXtract



## Biliary Extraction Balloon

### ScopeXtract - Extraction Balloon

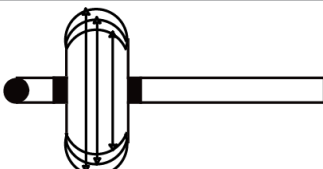	Reference	Ø Catheter (Fr / mm)	Length (cm)	Ø inflated balloon (mm)	Short guide wire compatible	QTY
	SEB-1012155451	7 - 5 / 2,3 - 1,7	200	10 - 12 - 15	✗	5
	SEB-1315185451	7 - 5 / 2,3 - 1,7	200	13 - 15 - 18	✗	5
	SEB-1012155453	7 - 5 / 2,3 - 1,7	200	10 - 12 - 15	✓	5
	SEB-1315185453	7 - 5 / 2,3 - 1,7	200	13 - 15 - 18	✓	5

## Characteristics

- ▶ The platoid shape of the balloon reduces the surface area in contact with the walls of the bile duct, so that:
  - to reduce friction and thus resistance to shrinkage,
  - to remove the majority of the stones in one pass.
- ▶ The balloon does not deform easily, thus avoiding lateral leakage during extraction.
- ▶ The conical design of the catheter tip allows for easy cannulation of the papilla.

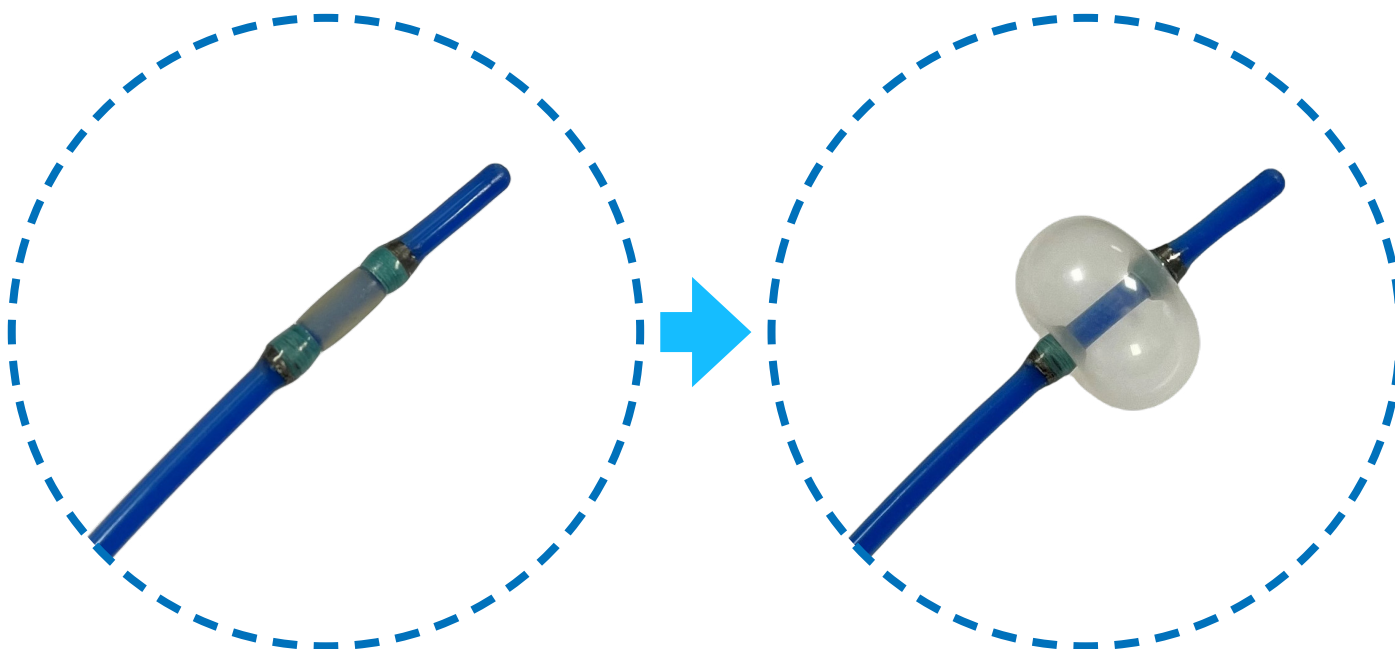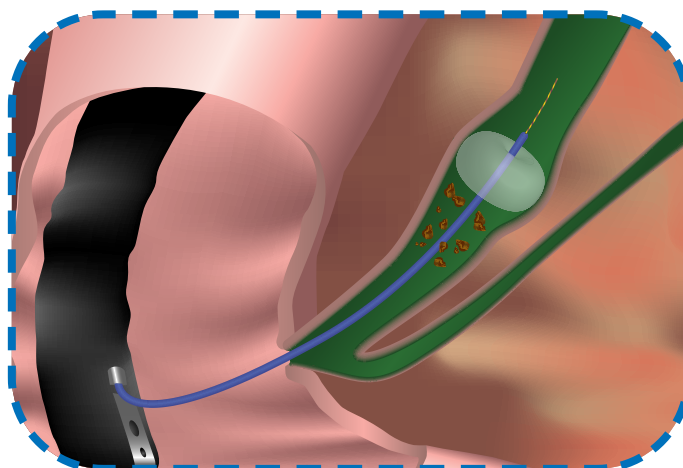

**Description:** ScopeXtract extraction balloon. **Indications:** Endoscopic removal of gallstones. **Class:** IIa. **Notified body:** CE0123. **Manufacturer:** Hangzhou AGS Medtech Co, Ltd. **Proper use:** Before use, carefully read the instructions in the package insert or on the device label. **Version:** SX\_E\_3. Rev.03. Date 15-07-2022

Creo Medical authorised local distributor

**Creo Medical UK Ltd**

Chestnut House, N.Minster Business Pk.

YO26 6QR York, United Kingdom

**Email:** info.uk@creomedical.com

**Tel:** +44 (0) 1904 786888

Technical changes and mistakes reserved. Images may show non standard options.

