

ScopeBasket

Biliary Extraction Basket

Designed for
endoscopic extraction
of gallstones and
foreign bodies

Key Characteristics

Nitinol model rotatable and with guidewire

4 nitinol single stranded wires ensure good expansion due to their shape memory.

The rotating spiral shape of the basket facilitates the grasping of gallstones.


Compatible with short guide wire and/or long guide wire.


Compatible with a universal lithotripsy handle in case of stone incarceration.

Standard model in stainless steel

4 single strands.


4 angled arms.

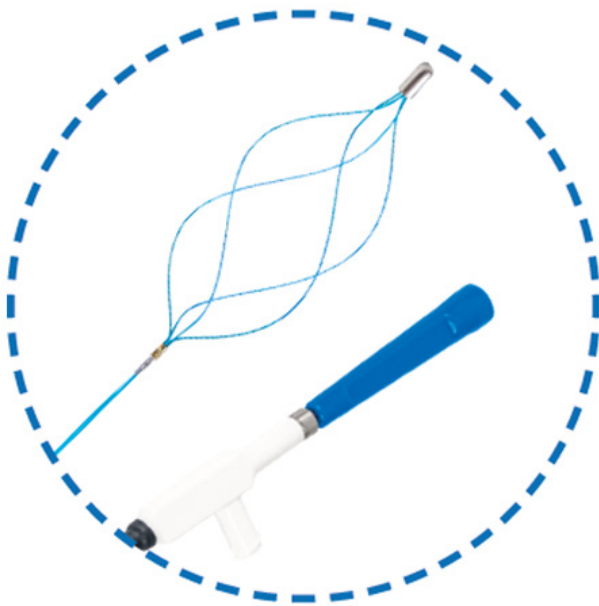

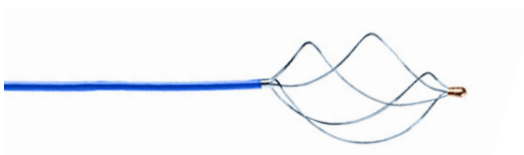

ScopeBasket - Rotatable (Nitinol) and wire-guided extraction basket

Reference	Ø Catheter (mm)	Length (cm)	Basket dimensions Width x Length	Short guide wire compatible	Recommended guide wire Ø	Box QTY
SEB-N2024G5062	3.2	195	20 x 40	Yes	0.025 / 0.035	5
SEB-N2524G5062	3.2	195	25 x 50	Yes	0.025 / 0.035	5
SEB-N3024G5062	3.2	195	30 x 60	Yes	0.025 / 0.035	5

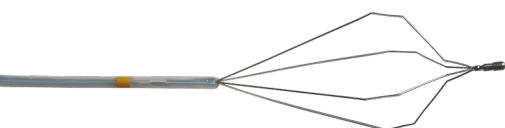

ScopeBasket - Standard (stainless steel) extraction basket

Reference	Ø Catheter (mm)	Length (cm)	Basket dimensions Width x Length	Short guide wire compatible	Recommended guide wire Ø	Box QTY
SEB-S2024G5063	3.2	195	20 x 24	No	0.025 / 0.035	5
SEB-S2524G5063	3.2	195	25 x 50	No	0.025 / 0.035	5

For orders and enquiries, please contact our Sales Team:

call: **+1 866-226-1170** email: **uscustomerservice@creomedical.com**

Creo Medical Inc

100 Reserve Road
Suite B400
Danbury
CT 06810

+1 866-226-1170

uscustomerservice@creomedical.com

in [showcase/creo-endoscopy](#)

in [company/creo-medical](#)

X [@CreoMedical](#)

X [@CreoEndoscopy](#)

Technical changes and mistakes reserved. Images may show non standard options.

Creo Medical, the CM logo, Kamaptive and Speedboat are trademarks of Creo Medical Limited. Certain names and marks are registered in the UK and other territories. All cited trademarks are the property of their respective owners. © 2026 Creo Medical. All rights reserved.

MKT-GEN-018-0 | 2026-October



**CREO
MEDICAL**